
# How to Use FlexPen (Insulin Pen)



### Step 1

#### Check the insulin

- a. Check label of pen and make sure it is the correct type of insulin
- b. Check for the expiration date
  - Never use insulin that is cloudy or not clear
- c. Remove pen cap (A)

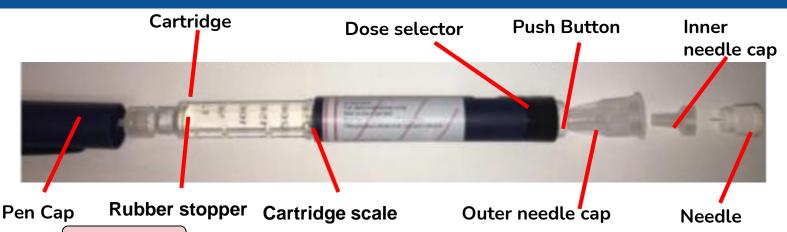

Α

## Step 2

#### Attach a new needle

- a. Remove the rubber stopper (B)
- b. Screw the needle tightly on
- c. Remove the outer needle cap (C)
- d. Remove the inner needle cap and dispose it (D)








## Step 3

## Priming the insulin needle

- a. Turn dose selector to 2 units (E)& Point needle straight upwards
- a. Check for air bubbles
- b. Press push-button all the way (F)
- c. Dose selector should return to 0 units
- d. A drop of insulin should appear at the needle tip

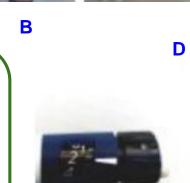





To learn more about our doctor's group, call Member Relations: (415) 590-7418



Ε

youtube.com/AAMGDoctors



# How to Use FlexPen (Insulin Pen)

### Step 4

#### Select the dose

 Adjust dose selector to the select the number of units you need to inject



G (28 units selected)

### Step 5

#### Administer the dose

- a. Clean injection site with alcohol pad and wait for it to dry
- b. Inject needle all the way into the skin (H)
- c. Press and hold push button to give dose (I)
- d. Leave needle under the skin for at least 6 seconds for insulin to be delivered before removing the needle





Н

# **Injection Sites**

- Rotate injection sites between injection spots
- Insulin is injected into the layer of fatty tissue under the skin, not the muscle (at least 1.5 inches away from the last injection)
- Do not inject into the same spot every time
- Do not inject into moles or scars



J

## Step 6

## After the injection

- Remove the needle and place into red sharps container (K)
- b. Replace the pen cap onto FlexPen
- c. Store FlexPen





To learn more about our doctor's group, call Member Relations: (415) 590-7418

#### aamgdoctors.com



youtube.com/AAMGDoctors facebook.com/DoctorsAAMG